

Mouse and Rat Intubation Kits



Features & Benefits:

- **Integrated anesthesia facemask**
Deliver isoflurane during intubation
- **Neck positioner**
Align airway structures
- **Fiber optic lighting**
Illuminate airway for enhanced visualization

Order your Intubation Kit:

888.572.8887



Kent Scientific
CORPORATION


The intubation of small animals continues to be a challenge for researchers. With an integrated anesthesia mask, the Kent Scientific intubation kit allows the delivery of isoflurane during intubation, making it easier to intubate small animals.

The ability to adjust the angle of the neck allows alignment of airway structures to visualize the vocal chords.

The fiber optic guide light illuminates the airway to enhance visualization.

Intubation kit includes:

- Intubation platform, positions animal at 60°
- Intubation stand with adjustable head positioning
- Anesthesia mask for anesthesia delivery during intubation
- Anesthesia mask tubing & connectors (36" length)
- Fiber optic light and charger
- Intubation catheters (5-pack)
- Suture, 1-meter length
- Small spatula to position animal's tongue
- Small mirror to confirm correct placement
- Power supply

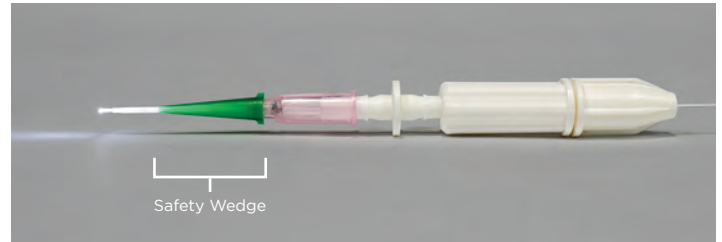

The Intubation Tube Safety Wedge is designed to assist in insertion of the catheter into the trachea and minimizes the chances of inserting the cannula too deeply. During intubation, the wedge will guide the catheter tip to the proper location between the larynx and carina. Once placed, the tight fit ensures bilateral ventilation and reduces air leaks from around the catheter.

1-Year 100% Parts & Labor Warranty



Integrate the SurgiSuite into your surgical set-up

Shown with optional FIR warming pad and controller, LED lighting and magnifier kit



Order your Intubation Kit:

888.572.8887

kentscientific.com



Kent Scientific
CORPORATION